



RAPID OPEN BITE CORRECTION (ROC)

TS^{*}₂

This innovative product has been designed by Dr. John Voudouris (University of Toronto and New York University). Open Bite are resolved in less than 4-6 months in conjunction with a light anterior box-pattern elastic. TS2 is placed palatally from upper 3-3 and/or lingually from lower 3-3 where it is totally invisible. The mildly pointed star-protrusions are comfortable reminders to efficiently prevent **Anterior Tongue Positioning**, while generous undercuts allow the use of intraoral elastics for the simultaneous crossbite correction. Completely invisible and esthetic lingual positioning that is generally comfortable.

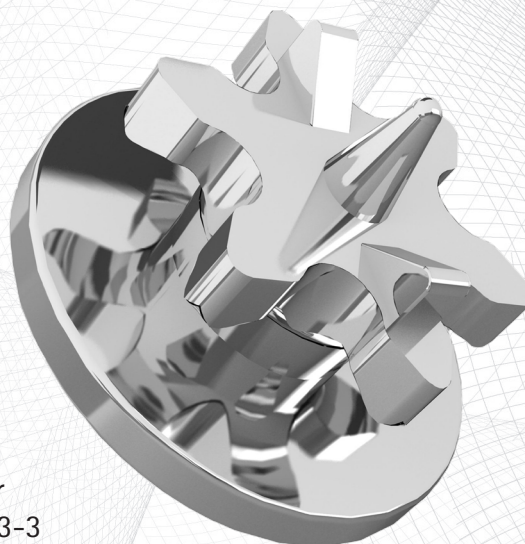


1. Nine rounded reminder protrusions prevents tongue thrusting (TS2 do not cut tongue, operator's glove or skin)

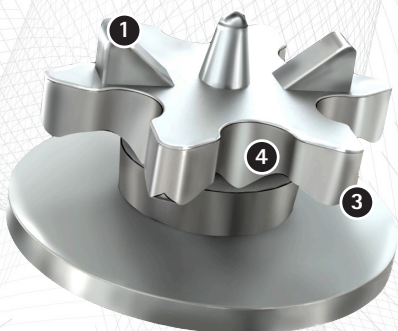

3. Six Undercuts for crossbite elastics in lateral open bite treatment that commonly require maxillary expansion

2. 80-gauge mesh bonding pad for increased shear bond strength

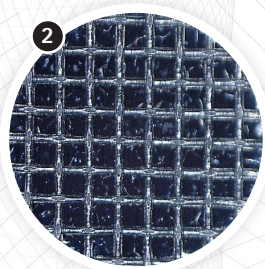

4. Mini tongue inhibitor button for good oral hygiene (alternative to Tongue crib) assisted by active self-ligation

CODE	PACKAGE
N30-TS2B	12 pcs

